

FlexiPump®


Dispensing pumps for liquid dispensing

Speed and precision in every drop




interscience

FlexiPump® is a range of peristaltic pumps for liquid dispensing. It adapts to all your needs, from media preparation to sample preparation. **By saving time, you can efficiently dispense all types of liquids**, even syrupy ones.

 **x 40**
in 1 min⁽¹⁾

 **x 29**
in 1 min⁽¹⁾

 **x 8**
in 1 min⁽¹⁾



Distribution accuracy and reproducibility

- **Accuracy** up to 0.5 %
- **Reproducible results** over long dosing distributions
- In compliance with ISO 6887 and ISO 7218 standards

PATENTED

Integrated scale ⁽²⁾

FlexiPump® Scale combines the power of a peristaltic pump with the precision of an weighing platform:

- **Integrated weighing platform** connected to the unit
- **Dose on scale** mode
- **Autocalibration with customizable frequency** during your dispensing series

Simplified distribution

- **Cross-contamination-free** liquid transfer
- **Compact** on the bench or under the PSM
- Fully autoclavable dispensing set with bottle
- Anti-drip system

Applications

- **Distribution:** Dilution tubes, bottles, jars, flasks, filter bags, sampling bags, ...
- **Pouring of Petri dishes**



Food



Veterinary



Pharmaceutical



Environmental



Public research



Cosmetics

⁽¹⁾ With **FlexiPump® Scale** in multi-dose mode with a delay of 0.1 s in optimal conditions.



FlexiPump®
Essential
functions



FlexiPump® Pro
High-performance
dispensing



FlexiPump® scale
Integrated
weighing platform

Flexibility and adaptability

The pumps that fits all your dispensing needs



4 dispensing modes

Continuous
Dose
Multi-dose
Dose on scale⁽²⁾





Wide range of volumes

Dispense from 50 µL to 99 L
in all types of containers
(plates, tubes, bottles,
filter bags, sampling bags)



3 dispensing options

Handy gun - ref. 561 001
Adjustable dispensing arm
ref. 561 002
Foot pedal - ref. 507 008

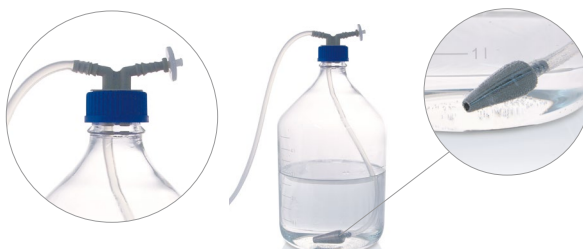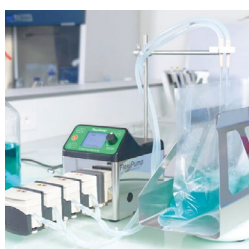

Jumbo Pack

Flow rates of up to 5L/min

BagOpen 3500 (inox) - ref. 563 350

Adjustable distribution system
ref. 561 002

Dispensing assembly + Y
(GL45 - Ø 6,4 mm) - 2x ref. 562 104

Booster kit - 2x ref. 511 000



**Tubings
guide**

Wide range of tubings

From Ø 1.6 to 8 mm to adapt
flow/accuracy to all your
applications.



Traceability

All distribution parameters
are traceable.

Printer⁽³⁾ - ref. 410 150

⁽²⁾ Only available on **FlexiPump® Scale** ⁽³⁾ Not available on **FlexiPump®**

Technical specifications

	FlexiPump®	FlexiPump® Pro	FlexiPump® scale
Reference	561 000	562 000	563 000
OVERVIEW			
All stainless steel 304 L	✓	✓	✓
Handy gun compatible	✓	✓	✓
Traceability (printer, PC, USB stick)	-	✓	✓
Double-head pump for high accuracy and high flow	-	✓	✓
Integrated weighing platform	-	-	✓
Autocalibration with customizable frequency (PATENTED)	-	-	✓
Jumbo pack compatible	-	-	✓
SPECIFICATIONS			
Volume dispensed in dose and multi-dose mode	50 µL to 99 L		
Volume dispensed in dose on scale mode	-		1 mL to 6.5 L
Weighing range	-		0.1 to 6500 g
Weighing resolution up to 200 g	-		0.01 g
Weighing accuracy	-		0 to 100 g: ± 0.05 g / 100 to 6500 g: 0.1 %
Typical dispensing accuracy with inner Ø 3.2 mm tubing	9 mL: ± 1.8% - 18 mL: ± 1% - 225 mL: ± 0.5%	9 mL: ± 0.9% - 18 mL: ± 0.5% - 225 mL: ± 0.4%	
Accuracy range in dose on scale mode	-		0 to 1% in accurate mode 0 to 2% in fast and normal modes (measured for each dispensing)
Dispensing modes	Continuous - Dose - Multi-dose		Continuous - Dose with/without autocalibration - Multi-dose with/without autocalibration - Dose on scale
Number of customizable programs	32		
Multi-dose pause	Adjustable from 0.1 to 99.9 s		
Watson Marlow™ pump	Single pump head	Low pulsation double pump head	
Flow rates	1.5 mL/min to 1.25 L/min	3 mL/min to 2.5 L/min	6 mL/min to 5 L/min ⁽⁴⁾
Pump head rotating speed	10 to 550 rpm		
Tubing material	Tygon™ (for pump head) and silicon		
Connections	Jack 3.5 mm	2 x USB type A, USB type B, 3.5 mm jack	2 x USB type A, USB type B, 3.5 mm jack, weighing platform socket
Voltage - Frequencies	100-240 V~ 50-60 Hz		
Warranty	1 year		
Spare parts availability	10 years		
In compliance with	ISO 7218, ISO 6887 and FDA BAM (Bacteriological Analytical Manual)		
TRACEABILITY			
Traceability	-	Operator name, dispensing mode accuracy, last calibration, date, time, speed, Ø tubing	Operator name, target weight, measured weight, dispensing mode accuracy, last calibration, date, time, speed, Ø tubing
Links	-	USB type A (printer, USB stick) USB type B (PC)	USB type A (printer, USB stick) USB type B (PC)
Data export	-	Compatible CSV Excel™, Open Office Calc™, LIMS	
WEIGHT AND DIMENSIONS			
Dimensions (w x d x h)	16.1 x 27.3 x 20.8 cm	16.1 x 33.2 x 20.8 cm	16.1 x 33.2 x 20.8 cm + Weighing platform: Ø 17 cm x 5.4 cm
Net weight	5.40 kg	5.55 kg	7.50 kg
Box dimensions (w x d x h)	50 x 40 x 40 cm		
Gross weight	8.20 kg	8.45 kg	11.26 kg

Delivered with 1 or 2⁽³⁾ Watson Marlow™ peristaltic pump heads, 1 weighing platform, 1 x 3.2 mm dispensing set (nozzle + GL45 screw cap + filter + tubing + tube weight), 1 x 6.4 mm dispensing set⁽³⁾, 1 user manual, 1 power cord, 1 USB cable⁽³⁾, 1 free pack of blender bags, 1 certificate of conformity / warranty card.

(3) not available on FlexiPump® (4) with Jumbo pack

Certified production



Product made for INTERSCIENCE by Interlab, an ISO 9001 certified company.

Your local distributor

interscience

PARIS

Phone: +33 (0)1 34 62 62 61 - Email: info@interscience.com

FRANKFURT

Phone: +49 611 7238 7770 - Email: sales.germany@interscience.com

BOSTON

Phone: +1 781 937 0007 - Email: sales.usa@intersciencelab.com

SHANGHAI

Phone: +86 (0)21-64739390 - Email: sales.china@interscience.cn

SINGAPORE

Phone: +65 6977 7232 - Email: sales.asia@interscience.com

TOKYO

Phone: +81 3 6712 9715 - Email: sales.japan@interscience.com