


BagMixer® 400 S

400 mL lab blender

Ref 025 000



Standards:

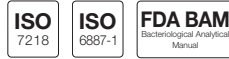

Manufactured under:



Product made for INTERSCIENCE by Interlab, an ISO 9001 certified company

BagMixer® 400 S is the next generation of lab blender: silent, easy-to-use and robust. With great new features, along with the silent blending, it is adapted to all kinds of applications and with a guarantee of optimal bacterial extraction, it is a great tool for lab analyses avoiding all risks of cross-contamination.

TECHNICAL SPECIFICATIONS

- All stainless steel 304L outer shell
 - Silent blending: less than 48 db
 - Noise cancelling casing
 - Multi-function digital screen
 - Useful volume: 50 - 400 mL
 - Variable speed (4 - 6 - 8 - 10 strokes/s)
 - Variable blending time (1 s - 1 h or ∞)
 - Adjustable blending power (patented)
 - Click & Clean® removable paddles (patented)
 - Side-by-side paddle stop
 - Q-Tight® closing: 250 kg pressure to avoid bag leaks from top
 - Lighted polycarbonate blending chamber with round corners: easy cleaning
 - Total Access system: total opening of the door for full access of the chamber (patented)
 - Super silent brushless motor
 - Integrated security drip tray
 - Paulstra™ feet: great stability during blending / non-marking Paulstra™ feet
 - Blending Light Code: green for Ready, orange for blending, red for leak detected
 - EnergySaver: 50% less energy consumed with brushless motor
 - Auto Sleep function: switch to sleeping mode after 30 mn of inactivity
 - Ambidextrous handle
-
- Voltage and frequency: 100-240V~ 50-60Hz
 - Fuses: T 3.15A L 240V 5x20 mm
 - Max power: 55 W
-
- Dimensions (w x d x h): 27.5 x 46.5 x 36 cm, weight: 25.2 kg
 - Box (w x d x h): 56 x 41 x 50 cm, weight: 28.6 kg
-
- 3-year warranty
 - Shock absorbers life-time warranty
-
- In compliance with: ISO 7218, ISO 6887-1 and FDA BAM (Bacteriological Analytical Manual)
 - Manufactured under CE, RoHS, WEEE - Products made for INTERSCIENCE by Interlab, an ISO 9001 certified company
 - Designed and made in France

DELIVERED WITH

- 1 integrated security drip tray
- 1 power cable
- 1 user's manual
- 1 quick user guide
- 1 free pack of blender bags