


dataLink® pro

Traceability system with semi-automated labelling

Ref. 439050

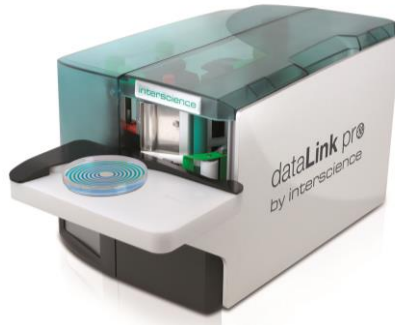

Certification :
in compliance with EU regulations

Traceability and High Technology:



Manufacturing:



dataLink® pro is a traceability system with semi-automated labelling for seeded Petri dishes. The operator prints the information of his choice on a label (volume, dilution, time, deposit mode...) which is stuck on the edge of the dish to be identified. Read in a **Scan®** or **ScanStation®** counter, these data allow the automatic configuration of the incubation and/or the counting of the colonies on the corresponding Petri dish.

TECHNICAL SPECIFICATIONS

- Thermal transfer printing
- LCD screen
- Movable cover
- **Semi-automated labeling tray**
- Synthetic rubber printing roller

- **Print resolution: 600 dpi**
- Max. printing speed: 150 mm/s
- Max. printing width: 105.7 mm
- Thermal transfer ribbon
- Max. roll diameter: 90 mm
- Max. length: 600 mm
- Width: from 25 to 114 mm

- **Pre-peel labelling**
- Centered material guide
- Mandrel diameter: 25 mm
- Label detection cell with position indicator
- **Label size: 5 x 50 mm, thickness: 0.03 - 0.6 mm**
- **Managing software: easySpiral®**

- Processor 32 bits, frequency: 800 MHz
- Random Access Memory (RAM): 256 MB
- Internal Flash Memory IFFS: 50 MB
- SD card slot (SDHC, SDXC): up to 512 GB

- USB 2.0 Hi-Speed port: computer connection
- Front USB port
- 2 USB ports
- SD reader
- Ethernet port
- RS232 serial port
- I/O port

- Power supply voltage and frequency: 100-240V ~ 50-60Hz
- Power: standby < 10W / nominal 100W
- Dimensions (W x D x H): 25.2 x 28.8 x 46 cm, net weight: 12 kg
- Packaging (W x D x H): 70 x 35 x 40 cm, gross weight: 16 kg

- Warranty: 1 year, after registration of the warranty card
- Made in Germany
- Complies with EU regulations
- **Compatible with interscience easySpiral®, Scan® and ScanStation® devices**

DELIVERED WITH

- 1 power cord
- 1 USB cable
- 2 user's manuals
- 1 thermal transfer ribbon
- 1 mandrel
- 2 rolls of 2500 adhesive labels
- 1 USB key to install software and drivers

COMPUTER SYSTEM REQUIREMENTS

- Processor: Intel i5, 2.8 GHz, quad core or higher
- Operating System: Windows™ 10 or Windows™ 11 (or higher)
- RAM: 4 GB minimum
- Hardware: 4 free USB ports
- Screen Resolution: 1280 x 1024 pixels or more

OPTIONS / ACCESSORIES

- Ref. 439 051 Lot of 3 rolls of 2500 adhesive thermal transfer labels: Printer stickers (5 x 50 mm) - Quantity: 3 x 2500
- Ref. 439 052 Thermal Transfer Film: Film for printer
- Ref. 522 000 Barcode reader