

Evaluation of the New Brilliance GBS Chromogenic Medium for Screening of *Streptococcus agalactiae* Vaginal Colonization in Pregnant Women

Paul O. Verhoeven,^{a,b} Pauline Noyel,^b Julie Bonneau,^b Anne Carricajo,^{a,b} Nathalie Fonsale,^b Alain Ros,^b Bruno Pozzetto,^{a,b} Florence Grattard^{a,b}

GIMAP EA 3064 (Groupe Immunité des Muqueuses et Agents Pathogènes), University of Lyon, Saint-Etienne, France^a; Laboratory of Bacteriology-Virology-Hygiene, University Hospital of Saint-Etienne, Saint-Etienne, France^b

Three commercial chromogenic agar media were evaluated for *Streptococcus agalactiae* screening in 200 vaginal swabs from pregnant women. The sensitivity and specificity were 94.3% and 100% for Granada medium (bioMérieux), 100% and 90.3% for Brilliance GBS medium (Thermo Fisher Scientific), and 100% and 98.8% for ChromID STRB medium (bioMérieux), respectively.

Streptococcus agalactiae, usually termed group B streptococcus (GBS), is one of the most important causes of early-onset neonatal infection (1, 2). The incidence of neonatal GBS infection ranges from 0.80 to 3.06 per 1,000 live births in developing countries (1). Guidelines recommend screening vaginal or rectovaginal GBS colonization in pregnant women at 35 to 37 weeks of gestation (3, 4). The prevalence of rectovaginal colonization varies from 6.5 to 36% in European countries (5). For pregnant women colonized by GBS, intrapartum administration of antibiotics is recommended to prevent GBS transmission to the newborn during delivery (3, 4).

The performance of microbiological methods for GBS screening has been greatly improved by the use of selective chromogenic media (6–9), together with rapid bacterial identification using matrix-assisted laser desorption ionization–time of flight mass spectrometry (MALDI-TOF MS) (10). The overnight broth pre-enrichment step, recommended in some countries (3, 4), increases the sensitivity of the GBS screening (8, 11), mostly in cases of low bacterial loads (11, 12). Enriched and direct cultures were shown to have similar performances (9), although failure to perform direct culture can decrease sensitivity due to overgrowth of competing organisms in enrichment broth (13). Finally, nucleic acid amplification tests (NAATs), notably those including the extraction step, reduce the delay in results and could be used at the admission of pregnant women in active labor in order to administer adequate antibiotic prophylaxis (14, 15). NAATs offer similar performances compared to direct selective culture (16), but controversial results were reported after comparison with enriched culture (7, 9, 17, 18). Recent guidelines advised the use of NAATs in selected circumstances, such as for term gestation women with unknown colonization status and no risk factors (3, 4).

The aim of this study was to assess the performance of the new Brilliance GBS chromogenic medium (Thermo Fisher Scientific, Dardilly, France) for the detection of *S. agalactiae* compared to those of two commercially available chromogenic media selective for this bacterium, ChromID STRB and Granada media (bioMérieux, Marcy l’Etoile, France).

From February to April 2013, 760 vaginal swabs (ESwab 480CE; Copan, Brescia, Italy) from pregnant women were sent to the microbiology department of the University Hospital of Saint-Etienne, France. The routine screening of GBS in pregnant

women was performed by direct culture following French recommendations (19). Approximately 50 μ l of ESwab medium was deposited on Granada medium by using the swab and spread with a sterile loop. Plates were incubated under anaerobic conditions at 36°C and read after 24 and 48 h of incubation following the manufacturer’s recommendations. Orange colonies were considered *S. agalactiae* without confirmation of identification. The prevalence of *S. agalactiae* vaginal colonization was 16.7% (127/760).

From this panel, 200 nonconsecutive samples were included prospectively in the study and plated onto Brilliance GBS, Granada, and ChromID STRB agar plates. Briefly, a 50- μ l volume of ESwab medium was plated onto each chromogenic plate by using the EasySpiral Dilute instrument (Interscience, Saint Nom la Bretèche, France) and samples were stored at 4°C. Plates were incubated at 36°C under aerobic conditions for GBS Brilliance and ChromID STRB media and under anaerobic conditions for Granada medium, as recommended by each manufacturer. All the plates were read after 22 to 24 h of incubation. Granada medium plates were read after 48 h in the routine workflow. Based on colony color, presumptive GBS colonies isolated on each chromogenic medium were systematically identified at the species level by MALDI-TOF MS (Microflex LT; Bruker) following the manufacturer’s recommendations. Non-target-colored colonies were also identified by MALDI-TOF MS if at least one chromogenic medium yielded GBS; otherwise, no MALDI-TOF MS identification was performed. Thirty-five (17.5%) of the 200 vaginal swabs were found positive for *S. agalactiae* by at least two chromogenic agar plates. No sample was found positive by only one medium. The performances of the chromogenic media are depicted in Table 1. The two GBS-positive samples missed by the Granada medium were found positive for nonhemolytic *S. agalactiae*. The 33 posi-

Received 19 October 2013 Returned for modification 16 November 2013

Accepted 2 January 2014

Published ahead of print 8 January 2014

Editor: W. M. Dunne, Jr.

Address correspondence to Bruno Pozzetto. Copyright © 2014, American Society for Microbiology. All Rights Reserved. doi:10.1128/JCM.02926-13

TABLE 1 Performances of Granada, Brilliance GBS, and ChromID STRB chromogenic media with 200 nonconsecutive vaginal swabs (including 35 positive ones by at least two culture methods) taken from women 35 to 37 weeks pregnant^e

Medium	No. of samples				Se (%)	Sp (%)	PPV (%)	NPV (%)
	Target color ^d		Nontarget color					
	GBS	Other sp.	GBS	Other sp.				
Granada	33	0	2 ^b	165	94.3	100	100	98.8
Brilliance GBS	35	16 ^c	0	149	100	90.3	68.6	100
ChromID STRB	35	2 ^d	0	163	100	98.8	94.6	100

^a Target-colored colonies were orange on Granada medium, pink on Brilliance GBS medium, and pink to red on ChromID STRB medium.

^b Growth of white colonies.

^c Six strains of *Streptococcus pneumoniae*, six strains of *Streptococcus mitis/oralis*, four strains of *Streptococcus salivarius*, and one strain of *Streptococcus parasanguinis*. One sample yielded pink colonies for both *S. pneumoniae* and *S. salivarius*.

^d One strain of *S. pneumoniae* and one strain of *Enterococcus faecium*.

^e Se, sensitivity; Sp, specificity; PPV, positive predictive value; NPV, negative predictive value; GBS, group B *Streptococcus*; sp, species.

tive samples on Granada medium were detected within 22 to 24 h of incubation; an extended incubation of 48 h failed to detect any further GBS-positive samples.

Additionally, the *S. agalactiae* load obtained on chromogenic media inoculated with the EasySpiral Dilute instrument was determined by automated reading using the Scan 1200 colony counter (Interscience) (Fig. 1A). A few samples yielding faintly colored

GBS colonies could not be reliably quantified from Brilliance GBS ($n = 3$), Granada ($n = 5$), and ChromID ($n = 4$) medium plates. After exclusion of these samples, the bacterial loads were available for the three media in 29 specimens, and these values ranged from 1.3 to 8.4 log CFU/ml with a strong correlation between media (Fig. 1B); the medians of the bacterial loads were 4.56, 4.82, and 4.57 log CFU/ml for Brilliance GBS, Granada, and ChromID

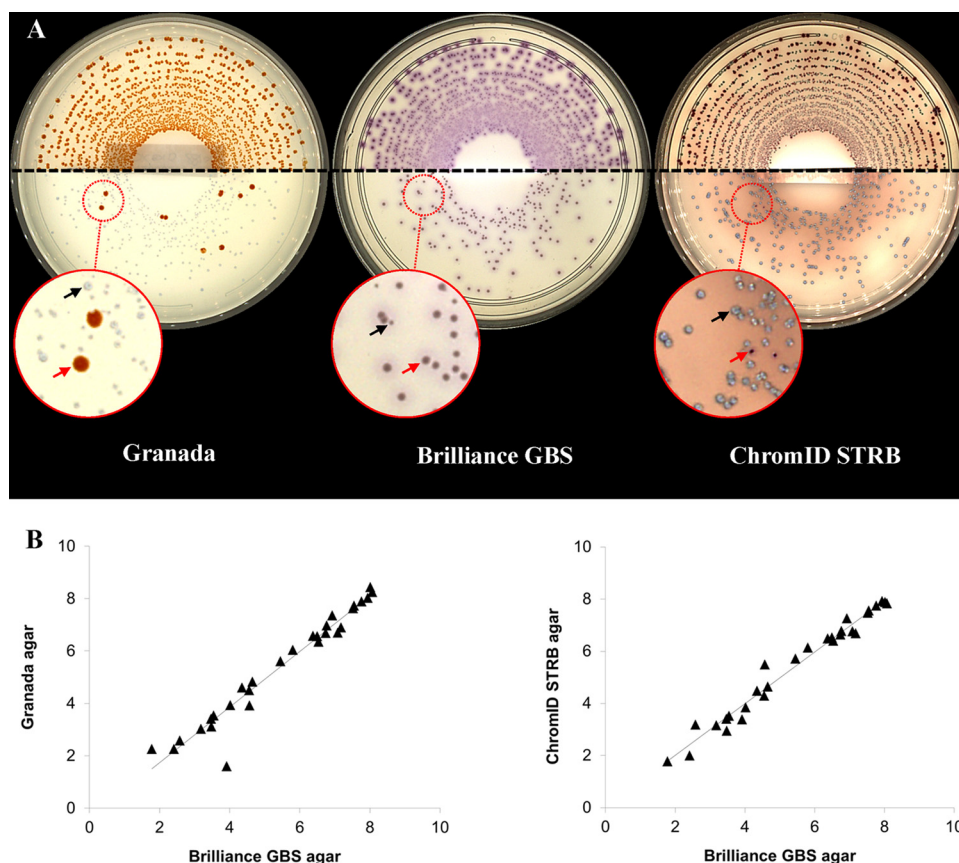


FIG 1 Comparison of three chromogenic media for GBS screening. (A) Pictures of Granada (left), Brilliance GBS (middle), and ChromID STRB (right) agar plates inoculated with two samples yielding averages of 6.42×10^7 CFU/ml (top) and 7.70×10^2 CFU/ml (bottom). In the enlarged parts of the picture, red arrows indicate target-colored colonies of *S. agalactiae* and black arrows non-GBS colonies. Agar plates were inoculated using the EasySpiral Dilute instrument and photographed with the Scan 1200 reader (see the text for details). (B) Correlations of GBS loads in 29 samples recovered on Granada and Brilliance GBS media (left graph) (Pearson coefficient of 0.97, $P < 0.0001$) and on ChromID STRB and Brilliance GBS media (right graph) (Pearson coefficient of 0.98, $P < 0.0001$).

STRB media, respectively (no significant difference by the Kruskal-Wallis test). Although it seems intuitive that a high level of GBS colonization in the vagina could favor neonatal infection, to our knowledge, the clinical significance of this parameter has not been evaluated prospectively. This evaluation could be achieved by using a chromogenic agar medium.

This study evaluated the new Brilliance GBS medium for screening GBS in vaginal samples of pregnant women on the basis of French recommendations. Although the three tested chromogenic media showed similar performances, the practicabilities differed greatly from one medium to another. The Granada medium failed to recognize 2 of the 35 GBS isolates that were shown to be nonhemolytic GBS (20). Plates were very easy to read because of the absence of atypical-colored colonies. The ChromID STRB medium showed excellent specificity and sensitivity but needed a high level of technical expertise due to the difficulty of recognizing the target colonies among all the colored colonies. The Brilliance GBS medium showed an excellent sensitivity and was easier to read than the ChromID STRB medium; however, it exhibited a lower specificity than the two other media because of the occurrence of false-positive results with other species of streptococci, highlighting the need to perform an additional identification step (agglutination test or MALDI-TOF MS) on colored colonies. Unlike Granada medium, Brilliance GBS and ChromID STRB media do not require an anaerobic atmosphere, which improves the laboratory workflow. From a clinical point of view, false-negative results could increase the risk of neonatal GBS infection due to the absence of antibiotic prophylaxis during labor. In contrast, false-positive results could lead to unnecessary administration of antibiotics during labor in a few women, which does not have a major adverse impact on mother or baby but leads to an overuse of antibiotics.

In conclusion, the new Brilliance GBS medium is a highly sensitive medium for GBS screening; it does not require anaerobic incubation. The positive colonies are easily recognized, but the identification of *S. agalactiae* at the species level would need an additional step to avoid false-positive results with other species of streptococci.

ACKNOWLEDGMENTS

This study was funded through a partnership with Thermo Fisher Scientific.

Emma Scopes is acknowledged for helpful discussion.

We declare having no conflict of interest related to this study.

REFERENCES

- Dagnew AF, Cunningham MC, Dube Q, Edwards MS, French N, Heyderman RS, Madhi SA, Slobod K, Clemens SAC. 2012. Variation in reported neonatal group B streptococcal disease incidence in developing countries. *Clin. Infect. Dis.* 55:91–102. <http://dx.doi.org/10.1093/cid/cis395>.
- Camacho-Gonzalez A, Spearman PW, Stoll BJ. 2013. Neonatal infectious diseases: evaluation of neonatal sepsis. *Pediatr. Clin. North Am.* 60:367–389. <http://dx.doi.org/10.1016/j.pcl.2012.12.003>.
- Verani JR, McGee L, Schrag SJ, Division of Bacterial Diseases National Center for Immunization and Respiratory Diseases, Centers for Disease Control and Prevention (CDC). 2010. Prevention of perinatal group B streptococcal disease—revised guidelines from CDC, 2010. *MMWR Recomm. Rep.* 59:1–36.
- Randis TM, Polin RA. 2012. Early-onset group B streptococcal sepsis: new recommendations from the Centres for Disease Control and Prevention. *Arch. Dis. Child. Fetal Neonatal Ed.* 97:F291–F294. <http://dx.doi.org/10.1136/archdischild-2011-300627>.
- Barcaite E, Bartusevicius A, Tameliene R, Klucinskas M, Maleckiene L, Nadisaukiene R. 2008. Prevalence of maternal group B streptococcal colonisation in European countries. *Acta Obstet. Gynecol. Scand.* 87:260–271. <http://dx.doi.org/10.1080/00016340801908759>.
- Tazi A, Réglie-Poupet H, Dautezac F, Raymond J, Poyart C. 2008. Comparative evaluation of Strepto B ID chromogenic medium and Granada media for the detection of group B streptococcus from vaginal samples of pregnant women. *J. Microbiol. Methods* 73:263–265. <http://dx.doi.org/10.1016/j.mimet.2008.02.024>.
- Smith D, Perry JD, Laine L, Galloway A, Gould FK. 2008. Comparison of BD GeneOhm real-time polymerase chain reaction with chromogenic and conventional culture methods for detection of group B Streptococcus in clinical samples. *Diagn. Microbiol. Infect. Dis.* 61:369–372. <http://dx.doi.org/10.1016/j.diagmicrobio.2008.03.007>.
- Craven RR, Weber CJ, Jennemann RA, Dunne WM, Jr. 2010. Evaluation of a chromogenic agar for detection of group B streptococcus in pregnant women. *J. Clin. Microbiol.* 48:3370–3371. <http://dx.doi.org/10.1128/JCM.00221-10>.
- El Aila NA, Tency I, Claeys G, Saerens B, Cools P, Verstraelen H, Temmerman M, Verhelst R, Vanechoutte M. 2010. Comparison of different sampling techniques and of different culture methods for detection of group B streptococcus carriage in pregnant women. *BMC Infect. Dis.* 10:285. <http://dx.doi.org/10.1186/1471-2334-10-285>.
- Cherkaoui A, Emonet S, Fernandez J, Schorderet D, Schrenzel J. 2011. Evaluation of matrix-assisted laser desorption/ionization-time of flight mass spectrometry for rapid identification of beta-hemolytic streptococci. *J. Clin. Microbiol.* 49:3004–3005. <http://dx.doi.org/10.1128/JCM.00240-11>.
- Louie L, Kotowich L, Meaney H, Vearncombe M, Simor AE. 2010. Evaluation of a new chromogenic medium (StrepB select) for detection of group B Streptococcus from vaginal-rectal specimens. *J. Clin. Microbiol.* 48:4602–4603. <http://dx.doi.org/10.1128/JCM.01168-10>.
- Jordan JA, Hall G, Davis T. 2010. Multicenter study evaluating performance of the Smart Group B Streptococcus (GBS) assay using an enrichment protocol for detecting GBS colonization in patients in the antepartum period. *J. Clin. Microbiol.* 48:3193–3197. <http://dx.doi.org/10.1128/JCM.00106-10>.
- Dunne WM, Jr, Holland-Staley CA. 1998. Comparison of NNA agar culture and selective broth culture for detection of group B streptococcal colonization in women. *J. Clin. Microbiol.* 36:2298–2300.
- Poncelet-Jasserand E, Forges F, Varlet M-N, Chaleur C, Seffert P, Siani C, Pozzetto B, Ros A. 2013. Reduction of the use of antimicrobial drugs following the rapid detection of Streptococcus agalactiae in the vagina at delivery by real-time PCR assay. *BJOG Int. J. Obstet. Gynaecol.* 120:1098–1108. <http://dx.doi.org/10.1111/1471-0528.12138>.
- El Helali N, Giovangrandi Y, Guyot K, Chevet K, Gutmann L, Durand-Zaleski I. 2012. Cost and effectiveness of intrapartum group B streptococcus polymerase chain reaction screening for term deliveries. *Obstet. Gynecol.* 119:822–829. <http://dx.doi.org/10.1097/AOG.0b013e31824b1461>.
- El Helali N, Nguyen J-C, Ly A, Giovangrandi Y, Trinquart L. 2009. Diagnostic accuracy of a rapid real-time polymerase chain reaction assay for universal intrapartum group B streptococcus screening. *Clin. Infect. Dis. Off. Publ. Infect. Dis. Soc. Am.* 49:417–423. <http://dx.doi.org/10.1086/600303>.
- Atkins KL, Atkinson RM, Shanks A, Parvin CA, Dunne WM, Gross G. 2006. Evaluation of polymerase chain reaction for group B streptococcus detection using an improved culture method. *Obstet. Gynecol.* 108:488–491. <http://dx.doi.org/10.1097/01.AOG.0000228961.42272.31>.
- Church DL, Baxter H, Lloyd T, Miller B, Gregson DB. 2011. Evaluation of the Xpert® group B streptococcus real-time polymerase chain reaction assay compared to StrepB Carrot Broth™ for the rapid intrapartum detection of group B streptococcus colonization. *Diagn. Microbiol. Infect. Dis.* 69:460–462. <http://dx.doi.org/10.1016/j.diagmicrobio.2010.11.002>.
- Haute Autorité de Santé (HAS). 2001. Antenatal prevention of the risk of early neonatal bacterial infection. Clinical practice guidelines. Haute Autorité de Santé, Saint-Denis, France. http://www.has-sante.fr/portail/upload/docs/application/pdf/Antenatal_prevention.pdf.
- Nickmans S, Verhoye E, Boel A, Van Vaerenbergh K, De Beenhouwer H. 2012. Possible solution to the problem of nonhemolytic group B streptococcus on Granada medium. *J. Clin. Microbiol.* 50:1132–1133. <http://dx.doi.org/10.1128/JCM.05372-11>.